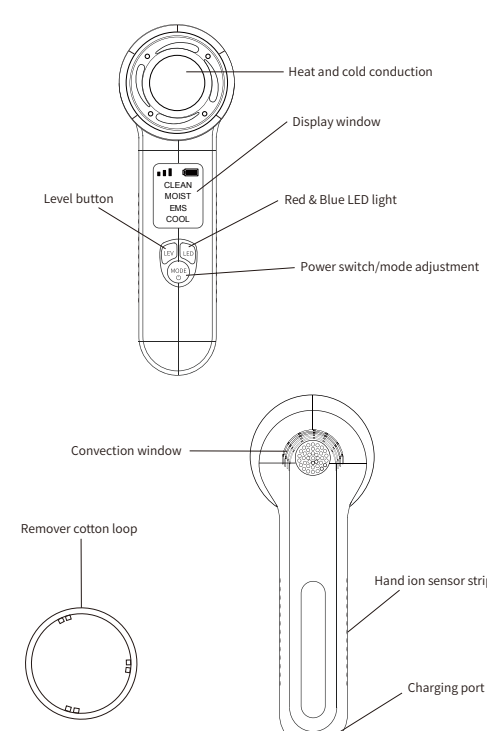


STRUCTURE COMPONENT	INTRODUCTION
 Heat and cold conduction Display window Level button Red & Blue LED light Power switch/mode adjustment Convection window Remover cotton loop Hand ion sensor strip Charging port	Warm effect - Warming, makes people feel warm and, helps to relax muscles, relieve pain; - Putting a sponge to deep cleansing, better able to expor the facetoins, blackheads, dirt. 3. EMS beauty principle (Electrical Muscle Stimulation) EMS is used to activate the skin cell, stimulating rejuvenation of collagen, to improve skinelasticity, soften wrinkles and leave skin looking smooth and healthy. 4. LED photon beauty principle - Red lightwavelength 620nm(± 10nm)18 LEDs deeply reaches into the 630nm skin to promote blood circulation and the collagen regeneration enhance the skin elasticity, to make the skin to be smooth and moisture. Be suitable for the slack & wrinkly skin. - Blue lightwavelength 470(± 10nm) 18LEDs deeply reaches into the 480nm skin, to kill the bacteria, promote the blood circulation, can penetrate into 0.25mm under the cortex, destroy the bacteria living environment, anti-inflammatory,bactericidal, remove acne, acne scars and so on effectively. 1.Acoustic vibrations beauty principle Each person's face has at least three million pores, while cleansing liquid and hands only clean skin surface, it is difficult to remove the oil pores,only clean out the pores of the oil residue, can the skin recover from the deep clean. Lessons learned 46 years specializes in skin care experience, the whole development process guided by skin specialists, combined with our advanced engineers to sonic, micro-vibration technology, fitting face of intelligent design, no dead comprehensive clean your skin.8000 times per minute, high-frequency vibration, vibrate out of deep pores oil and deeply clean every pores. Pores lightweight transparent, delicate skin smooth, once again change the skin more through clearlines. 2、Principles of hot and cold beauty A Japanese professor at the university shows that repeatedly reaching the skin at a temperature of 42℃ helps to produce HSP protein, repair skin cells and strengthen skin elasticity to resist the formation of wrinkles. The warm 42℃ treatment helps to dilute the pores for deep cleansing, while the cold treatment with a cold effect of up to 6℃ instantly tightens pores and locks in nutrition. Alternating warm and frozen skin care, it stimulates the liquid circulation at the bottom of the skin and speeds up the metabolism of skin cells. Besides effectively reducing the elimination of eye circles and single eye lines, it also brighten skin. Even without makeup, the skin is also revealed. Combined with micro-vibration, it is beneficial to enhance the skin's firmness and relieve fatigue. Ice sense effect - Help skin firmness, maintain flexibility,prevent allergies and other problems - Shrink pores, lock moisture and wrinkle. - Can soothe muscles, relieve skin fatigue.