




Multi head fascia gun massager manual

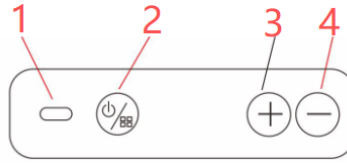

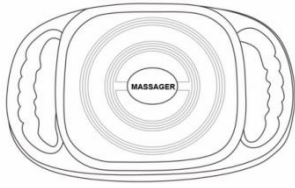

Model: KM04

Operation panel and button instructions:

1. Type-C charging port, 2. On/Off button, mode switch button, 3. Increased power, 4. Reduced power

 Button: Long press the button to turn on/off the massager. In power on mode, short press the button to enter mode selection, with 1-9 modes to choose from. Short press once each time to cycle between different modes, and display the usage time during use.

 Key: In the power on state, press the key to increase the force,  Key: In the power on state, press the key to reduce the force, display L1-L9, and return to the machine countdown in 2 seconds.

When the battery is low, you can use a universal 5V 2A charger to charge this product. The power indicator light flashes white during charging, and the prompt light stays on after being fully charged. The charging time is about 5 hours. When not using rechargeable batteries for a long time, it is recommended to charge them at least once a month. When the charging plug is inserted into the charging port, the machine cannot operate. It can only be used after unplugging the charging plug.

Common problem handling methods:

If you encounter any problems while using this product, please refer to the following tips and possible solutions.

Users are not allowed to disassemble or repair this product by themselves. Touching internal components may cause malfunctions or electric shock accidents.

1. The machine is not working: ● Please confirm that the switch is turned on. ● Please check if the battery is charged.
2. When working, suddenly stop: ● You may accidentally press the button on the controller. ● The 10 minute time set by the automatic timing device has passed. ● If the overheat protection device is activated, please turn off the main power of the product and wait for the massager to cool down before using it.

Precautions:

- If the machine fails to work or operates abnormally, it should be stopped immediately.
- Do not use for a long time, do not exceed 30 minutes.
- The product has a heating surface, and people who are not sensitive to heat must pay attention when using it.
- When unplugging or plugging the charging cable, do not forcefully pull or tug on the connecting cable.
- Do not use in damp environments, and do not place this product in water or other liquids.
- Before discarding this product, follow the steps below to remove the battery: cut off the power supply of the product, find the location where the battery is placed, and then use scissors to cut off the battery leads one by one before removing the battery. Dispose of the battery safely

Warning:

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

 Made in China