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40KHz Cavitation Slimming Machine

User Manual

QYX-CT20

持有认证：CE, FCC, PSE

Certifications: CE, FCC, PSE

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1. Overview

The 40KHz Cavitation is the ultrasonic instrument which gathered the strong sound wave head strong acoustic emission frequency of 40000Hz sound waves, after entering the body can make body fat cells to produce a strong impact and friction of motion between fat cells, can heat consumption effectively, water consumption in cells, and then fat cells was reduced, and when the sound wave vibration can make the cell produces strong crack, instantaneous blast cells, fat cells decreased, so as to remove fat effect.

2. Product features

- ★Cavitation Frequency: 40KHz -- Low Frequency.
(Strong penetration of ultrasonic)
- ★High-Intelligent Touch Screen Operation.
- ★Suitcase designed. It's portable and easy to carry.
- ★Time adjustment / intensity adjustment.

3. Specifications

Screen	8"LED touch color screen
Cavitation Head	X 1
Input Voltage	AC 100V~240V, 50Hz/60Hz
Working Voltage	DC 12V
Consumption Power	10W~33W
Cavitation Frequency	40KHz
Power	10~33W
Dimension	360(L)*240(W)*105(H)mm
G.W	4 KGS

4. Operation panel pictures

Switch

Turn on the power switch, start up the device

Press into the next page

Please choose Chinese/English in the bottom right, then press the pink/blue button into the next page

Main operation page

The middle shows:

Intensity control increase/decrease touch keys.

Time control increase/decrease touch keys.

The upper left shows:

The internal vibration frequency setting, please don't change it wilfully (Upper: baseboard frequency -- fixed 13Hz; Lower: vibration frequency -- depends on our factory setting).

The upper right shows:

The turn key -- turn back to the language selection page.

The lower right shows:

Start/stop touch keys.

5. Introduction of operation steps

(1) Connect power supply, turn on the power switch, the device start-up.

(2) LCD display language choice.

Press the left / right button to choose in English / Chinese language, and then enter the next page after confirmation.

(3) After entering the main operation panel, press the function button selection operation as follow:

① The middle shows:

Intensity control increase/decrease touch button.

Vibration amplitude is divided into 9 levels.

According to the customer bear ability to adjust, we suggest to adjust 6-8 levels.

Time control increase/decrease touch keys.

We suggest to adjust :3-5 minutes every in early use.

② The upper left shows:

Do not change the instrument internal vibrations frequency setting.

Upper: baseboard frequency fixed: 13Hz.

Lower: vibration frequency, depends on our factory setting.

③ The upper right shows:


The return key, turn back to the language selection screen.

④ The lower right shows: start / stop touch keys.

6. Use position: Body


(1) The thigh

- ① Inner thigh→thigh→outer thigh pull upward to the groin, as figure A, B shows.
- ② Touch inner and outer thigh back and forth in Z or eight methods as figure C, D shows.

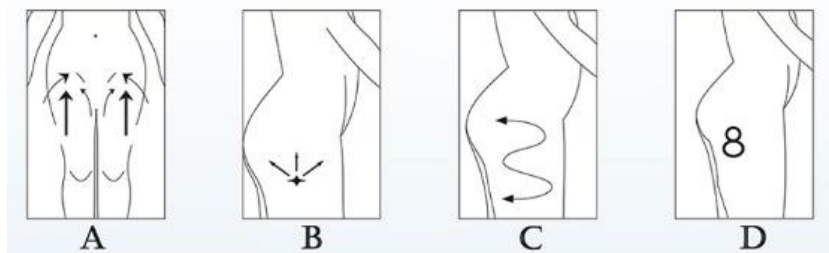(2) Abdomen

- ① Use from the waist move to the middle in circle way, as figure E shows.
- ② Use from the waist and the groin back and forth, as figure F shows.

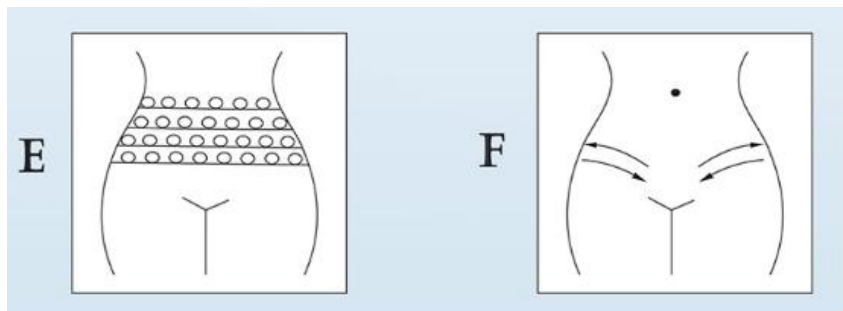7. Caution

- (1) High precision instrument which avoids placing in the wet location.
- (2) Instrument can't be hit
- (3) It is prohibited from using as below:
 - ① The wire is broken or getting hot.
 - ② Socket is loose or shake.
 - ③ It is forbidden to damage, processing, bend and twist plug too much.
 - ④ Do not pull the plug while operating.
 - ⑤ It is prohibited to use other power cord at the same time.
 - ⑥ Keep out of reach of children/ infant.
 - ⑦ It is prohibited to use chemical to clean the instruments.
 - ⑧ Do not pull out the plug with wet hands.
- (4) People who can not use
 - ① Serious diseases such as heart disease, gastric ulcer, duodenal ulcer, serious tummy bug, kidney is bad, etc.
 - ② It is strictly prohibited for lactating women, abdominal surgery under three months or menstrual period .
- (5) This instrument is limited to use around the abdomen and the thigh. It is forbidden to use in the bone joints and face, etc.
- (6) If the instrument set 7-9 levels(continuous), we suggest that use time no more than 3-5 minutes (in accordance with the individual acceptance), because the higher the vibration, the higher the temperature.
- (7) Contact with the skin directly, use personal skin lubricating

oil.

8. Maintenance

- (1) Please note that it is important to pay attention to maintenance, so as not to cause the shorten product life.
- (2) Cavitation head should be handled with great care, it is noted that using the device with hold the handle, do not pull the output line by intention (it can cause loose part inside the interface, causing poor contact and short circuit).
- (3) If it works 30minutes continuously, after that, we suggest that the head should be stopped in case of overheat.

9. Warranty

CE approval. No.: 14KWE091916E

FCC approval. No.: 14KWE091919F

PSE approval. No.: 14KWE091920J

- (1) Warranty: 12 months.
- (2) For the following user personal causes of the fault, does not provide free warranty service:
 - ① Disassembling or modified the device without authorization.
 - ② The external force damage, or natural disasters, human factors damage.
- (3) No prior notice will be given for any change of products. Please do not hesitate to contact us if you have any further question.

Maintenance		
Product:	Purchasing date:	Expiration date:
Distributor:	Address:	Cell phone:
Client:	Address:	Cell phone:
Date	Fault cause	Content
Remarks:		